

MODEL

Jazzy 1450-2VHD-SS	Group 2 Very Heavy Duty: K0826
Jazzy 1450-2VHD-C	Group 2 Very Heavy Duty: K0827
Jazzy 1450-2SPVHD-SS	Single Power Very Heavy Duty: K0839

Weight Capacity	600 lbs.
Turning Radius	26.5"
Overall Width ²	29.25"
Overall Length ²	38.25"
Maximum Speed ¹	Up to 4.5 mph
Ground Clearance ²	1.5"
Front Wheels	4" solid anti-tips
Drive Wheels	14" solid
Rear Wheels	9" solid casters
Drivetrain	Two-motor, front-wheel drive
Braking System	Electronic, regenerative disc brakes
Suspension Type	Swing Arm
Standard Seat Size	24" x 24"
Per Charge Range ^{2,4}	Up to 15.7 miles
Standard Electronics	90A, PG VR2
Component Weights	Base: 195 lbs. (heaviest piece) Seat: 52 lbs. (high-back) Batteries: 52 lbs. each
Battery Requirements ^{3,5}	(2) 12 volt, Group 24
Battery Charger	8 amp, off-board
Warranty (K0826/K0827)	Frame: Lifetime limited Drive motors: 13-month Electronics: 13-month Batteries: 13-month
Warranty (K0839)	Frame: Lifetime limited Drive motors: 18-month Electronics: 2-year

¹ Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tire condition. These specifications can be subject to a variance of (+/- 10%).

² Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or - 3%).

³ AGM or gel cell type required.

⁴ Tested in accordance with ANSI/RESNA, WC Vol2, section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.

⁵ Battery weight may vary based on manufacturer.

NOTE: All specifications subject to change without notice. The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%. Special order seats or non-standard dimensions will not be fold-flat and depth-adjustable.

PRIDE[®]
MOBILITY
LIVE YOUR BEST[®]

(US) 800-800-8586 • (Canada) 888-570-1113

gojazzy.com |    

HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. DME providers are responsible for determining appropriate billing codes when submitting for insurance reimbursement.

The information contained herein is correct at time of publication; we reserve the right to alter specifications without prior notice.

0001706-071020

Jazzy[®]
1450

**ROBUST AND
HEAVY-DUTY
DESIGN**



AMERICA'S FAVORITE POWER CHAIR

gojazzy.com |    

HEAVY-DUTY PERFORMANCE



Depth-adjustable, limited-recline, fold-flat, high-back seat

Universal, fixed, in-line joystick mount

Oversized rear casters for optimal performance and comfort

Large 17" x 13" foot platform



Superb traction, flat-free, drive tires



Extremely efficient motor package


The **Jazzy® 1450** offers heavy-duty construction and performance components. This bariatric power chair features a 600 lb. weight capacity, front-wheel drive for superb indoor maneuverability in the tightest of spaces and large front drive wheels for extra absorption over obstacles. The Jazzy 1450 combines superb climbing capabilities with front anti-tips for transitioning over grade changes and handling various terrains easily.

- Depth-adjustable, limited-recline, fold-flat, high-back seat
- Flat-free drive tires
- Universal, fixed, in-line joystick mount
- Large 17" x 13" foot platform
- 600 lbs. weight capacity
- Up to 4.5 mph maximum speed

Available Colors



Viper
Blue

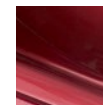

Candy
Apple Red

Accessories

- Lap belt
- Cane/crutch holder
- Cell phone holder
- Cup holder
- Flag
- Oxygen tank holder
- Rear basket
- Saddle bag
- Walker holder
- Weather cover
- XLR USB charger
- Swing-away joystick
- Military patches

Options

- 17"x16" foot platform
- HD drop-in legrests
- HD elevating legrests
- Swing-away joystick

Optional Seating

- High-back, black vinyl (pinchless hinge)
Available sizes:
22" W x 20" D
24" W x 20" D
24" W x 22" D
24" W x 24" D
- High-back, solid seat pan (pinchless hinge)
Available sizes:
22" W x 20" D
22" W x 22" D
- Heavy-duty high-back, black vinyl with headrest
Available sizes:
26" W x 20" D
26" W x 22" D